Applying Ultrasound Gel to the Mouthpiece(s)

The Aevo Oral Ultrasound Gel is intended to be used with the Aevo System™ device. Please refer to the Aevo System Instructions for Use provided with the Aevo System™ device or on the Aevo System™ website.

 www.aevosystem.com/instructions-for-use/

WARNING: Before use, ensure that the ultrasound gel is not expired and the gel pouch is not damaged. Only use the ultrasound gel provided with the Aevo System™. Other ultrasound coupling gels available on the market are not designed or tested for use with the Aevo System™.

Serious incidents are reported to the manufacturer and the competent authority. In case of serious incident that has happened in relation to the Aevo Gel, please report to SmileSonica Inc.

Before starting a treatment, tear open a corner of the gel pouch. All the gel from the pouch should be applied as a thin layer (3–4 mm thick) to the inside of the mouthpiece (see figure). This ensures contact between the mouthpiece and your gums.


Do not swallow the gel. Rinse out your mouth with water after each treatment.

NOTE: Use the ultrasound gel during every treatment or the treatment may not be effective.

IMPORTANT:

- Gel pouches are single use. This prevents contamination or infection.
- Store gel at room temperature between 15°C and 30°C.
- Expiry date format on gel pouches is YYYY-MM-DD.

Aevo System

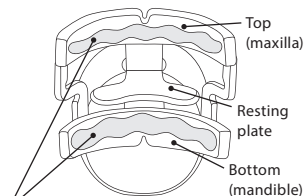Apply gel to the mouthpiece sidewalls

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**Emergo Europe B.V.**